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Mi, **Momert Zrt.**

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## OPERATING MANUAL ELECTRONIC BABY AND CHILD SCALE MODEL 6475

Thank you for choosing our product. Your electronic baby scale is a sensitive measuring instrument, therefore please protect it from overload and from big dynamic loads! Never leave your baby unattended on the tray!

Please, read the instruction for use before using the scale. Your new electronic baby scale has been designed for measuring your baby's weight and displaying his/her weight difference.

### TECHNICAL SPECIFICATION:

- Automatic hold function
- Precision measuring system
- Weight difference display function
- Tare function
- Low battery indicator ("Lo")
- Auto power off

Technical parameters	
<b>Voltage</b>	2 pcs 3V CR2032 lithium batteries
<b>Weighing range</b>	20 kg
<b>Weighing increments</b>	10 g
<b>Weighing units</b>	Kg, Lb, St
<b>Working/Storage temperature</b>	20-25°C/15-40°C
<b>Display size</b>	LCD diameter= 8,0 cm , 3,15 inch


### IMPORTANT SAFETY INSTRUCTIONS:

- Make sure you follow the marked polarity of battery when inserting and replacing it.
- If the appliance will be out of use for a longer time, remove the battery.
- Do not leave the unit unattended when in use.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
- Use the product only for the intended purpose specified in the Operating Manual.
- Use the unit only on a stable horizontal surface.
- Before cleaning, assembly, or disassembly of the cover, switch the unit off and remove the battery!
- Do not use abrasive or chemically aggressive substances for cleaning the unit.
- Do not use the unit if it does not operate properly, or if it has been dropped, damaged or exposed to any liquids. Take the unit to an authorized service center for testing and repair.
- The unit is suitable for domestic use only. It is not designed for commercial operation.
- Keep the unit away from heat sources, such as radiators, ovens etc. Do not expose the unit to direct sunlight or moisture.
- Do not touch the unit with wet or damp hands.
- The product is a fragile measuring device. Prevent it from crashing or falling.
- Do not repair the appliance yourself. Contact an authorized service center.

**FAILURE TO FOLLOW THE MANUFACTURER'S INSTRUCTIONS MAY LEAD TO REFUSAL OF WARRANTY REPAIR!**

### PRODUCT DESCRIPTION:

1. Weighing surface (baby tray)
2. ON/OFF button
3. TARE button
4. Kg/Lb/St button
5. Display
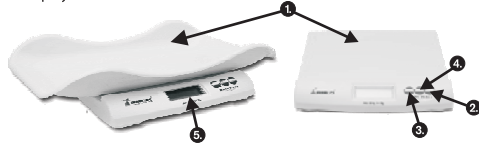
### OPERATING BUTTONS:

- The scale can be switched on and off with the ON/OFF button. (Push the button once.)
- With the TARE button you can use the tare function. "0.00" appears on the display 2 seconds after pressing the button. If you don't need the weight difference to be shown, then by pressing the TARE button for 3 seconds the memory can be deleted. "M" will disappear from the display.
- Kg, Lb, St button is suitable to choose the required unit of measurement, the actual unit of measurement can be seen on the right side of the display.

### USING THE SCALE:

1. These scales run on 2 pcs of replaceable lithium batteries. Remove the protective film before use (Fig. 1).
2. Before weighing a baby, snap the tray to the measuring surface so the tray would stand stable (Fig. 2).
3. There is no need for installing the tray if you want to weight a child.
4. By pressing the ON/OFF button once "0.00" appears on the display (Fig. 3).
5. Place the baby or stand the child on the scale. The measured weight is going to be displayed in a few seconds, and stored in the memory automatically. An "M" is displayed indicating the presence of data in the memory. After the next weighing the screen displays the measured weight and the difference between the measured weight and the previous weight by turns of 2 seconds. The lastly measured weight is automatically stored in the memory in case of both manual and automatic switching off.

Fig. 1

Fig. 2

Fig. 3

### MEMORY FUNCTION:

After switching on the scale, an "M" is displayed which means that the last measured weight is stored in the memory, so that the weight difference could be displayed. If you don't need the weight difference to be shown, then by pressing the TARE button for 3 seconds the memory can be deleted. "M" will disappear from the display. If you take out the batteries, then the memory will also be deleted.

### ABOUT TARING:

If you use not only the tray and the added weight (blanket, pillow etc.) is lighter than 5kg, after pressing the TARE button the screen displays "0.00" and the maximal measurable weight will be 20kg. If the added weight is heavier than 5kg, after pressing the TARE button the screen displays "0.00" but the maximal measurable weight will be

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Ryc. 1

less than 20kg (more precisely: 20kg minus the added weight).

#### OTHER INFORMATION ON THE DISPLAY:

If the "Lo" warning is displayed after switching on (Fig. 4), or there is nothing on the display, the batteries are dead and has to be replaced.

Please open the batteries holder at the backside of the scale and carefully take out the batteries by pressing the small silver colour holding plate to get the batteries out of the batteries holder. Replace the low batteries to a new batteries. The low batteries should be thrown to the appropriate place.

#### BATTERIES ARE NOT INCLUDED IN WARRANTY!

If the sign "ERR" appears on the display (Fig. 5), it means that the scale is overloaded, in this case the weight must be reduced (an overload can cause sensor faults).

If the message "OUT2" appears on the display (Fig. 6), it means sensor fault (after overloading or after a sudden, big dynamic load). In this case, we have no warranty liability for the product.

#### CLEANING AND MAINTENANCE:

The surface of the unit may only be cleaned using a wet cloth. Do not use any detergents or hard objects as they may damage the unit's surface!

**Never clean the unit under running water, do not rinse it, and do not submerge it in water!**

Fig. 4

Fig. 5

Fig. 6

#### REPAIR AND MAINTENANCE:

Any extensive maintenance or repair that requires access to the inner parts of the product must be performed by an expert service centre.

#### ENVIRONMENTAL PROTECTION:

- Packaging materials and obsolete appliances should be recycled.
- The transport box may be disposed of as sorted waste.
- Polyethylene bags shall be handed over for recycling.

#### Appliance recycling at the end of its service life:

 A symbol on the product or its packaging indicates that this product should not go into household waste. It must be taken to the collection point of an electric and equipment recycling facility. By making sure this product is disposed of properly, you will help prevent the negative effects on the environment and human health that would otherwise result from inappropriate disposal of this product. You can learn more about recycling this product from your local authorities, a household waste disposal service or in the shop where you bought this product.

#### WARRANTY:

We give a guarantee for this product for 12 months. In case of defect of product, to get your reclamation inspected and solved, please contact the commercial unit, where the product was bought by you, with this instruction manual booklet (and the guarantee card included). You may also contact directly the manufacturer with your reclamation. In this case the expenses of sending back the product to the manufacturer will be charged to you.

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