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Mi, **Momert Zrt.**

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Declare that product:

Name: Momert, **Model:** 1760, **Description:** Infrared ear and forehead thermometer

Is in conformity with the following standards:

Electromagnetic compatibility directive
IEC/EN 60601-1 Medical electrical equipment-Part1: General requirements for safety.
IEC/EN 60601-1-2 Medical electrical equipment-Part2: Collateral standard: Electromagnetic compatibility-Requirements and tests.
MDD (93/42/EEC); EN12470-5:2003; ASTM E1965-98; EN980: 2003;
ISO 14971 The Classification according to; IIEC/EN 60601-1 sub-clauses 5;
Internally powered equipment; IPX0
ROHS 2002/95/EC

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Tamás Baritsa
Chief Commercial Officer

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GB
INFRARED EAR AND FOREHEAD
THERMOMETER
MODEL 1760

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Operation manual

INFRARED EAR and FOREHEAD THERMOMETER

Model 1760

Thank you for purchasing a Momert product. We hope you will be satisfied with our product through its service life. Please study the entire Operating Manual carefully before you start using the product. Keep the manual in a safe place for future reference. Make sure other people using the product are familiar with these instructions.

Introduction

Utilizing infrared technology, this thermometer takes temperatures in seconds by measuring heat generated by measuring heat generated by the ear canal, the surface skin of the forehead or other objects. This product conforms to the provisions of the EC directive MDD(93/42/EEC). Its advantages include:

1. 4 in 1 functions: Ear/Forehead/Ambient temperature/Clock time
2. Probe Cover Free
3. Waterproof Probe Design
4. One-second Reading
5. Auto power off for power saving
6. Low-battery Indicator
7. Warning indication Indications for battery condition and measuring range.
8. Large LCD Display

Important Information Before Use

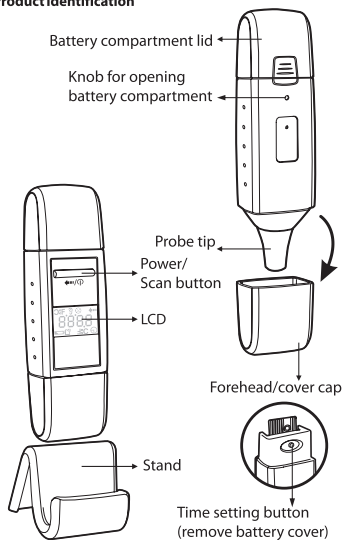
When using this product, please be sure to follow all the notes listed below. Any action against these notices may cause injury or affect the accuracy.

1. Do not disassemble, repair, or remodel the thermometer.
2. Be sure to clean the thermometer lens each time after usage.
3. Avoid direct finger contact with the lens.
4. It is recommended that user may take 3 temperatures. If they are different, use the highest reading.
5. Do not expose the thermometer to extreme temperature, very high humidity, or direct sunlight.
6. Avoid extreme shock or dropping the device.
7. Before the measurement, patients and thermometer should stay in steady state room condition for at least 30 minutes.
8. Avoid measuring temperature in 30 minutes after exercise, bathing, or returning from outdoor.
9. To protect the environment, dispose of empty batteries at appropriate collection sites according to national or local regulations.
10. Please use the thermometer solely for its intended purpose.
11. Allow one minute between successive measurements as slight variations may occur if measurements are taken over a short period of time.
12. There are no absolute body temperature standards. Keep reliable records of your personal temperature to serve as a reference for judging a fever.
13. Under any circumstances, the temperature taking result is ONLY for reference. Before taking any medical action, please

consult your physician.

14. If the appliance will be out of use for a longer time, remove the battery.
15. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
16. Use the product only for the intended purpose specified in the Operating Manual.
17. Before cleaning, assembly, or disassembly of the cover, switch the unit off and remove the battery!

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
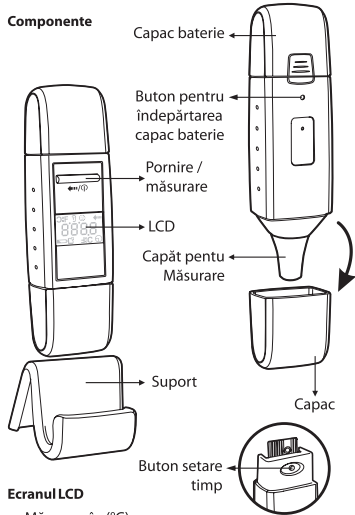
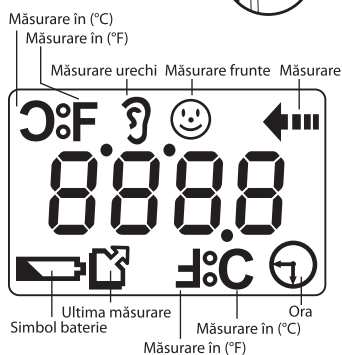
Description of LCD Display

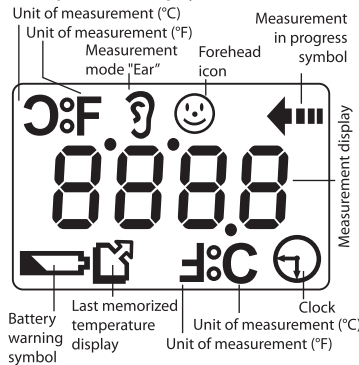
Time and Ambient Temperature

When the thermometer is used for the first time or after a battery change, the time must first be set. Then the ambient temperature measuring mode is automatically activated. If you do not need the room temperature or time to be displayed, or if you wish to preserve the battery, do not program in the time after changing the battery. This way you can turn on the thermometer for ear/forehead measurement as described. However, it then switches off completely.

Setting the time

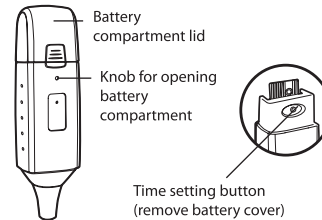
1. Open the battery compartment by pressing the knob for opening the battery compartment with a pointed object and removing the battery compartment cover. Press the time setting button for setting the time until the hour display flashes. The hours are set by pressing the button several times until the required hour setting is reached.

2. You can now switch from the hour setting to the minute setting by pressing the time setting button again. The minute display flashes. Adjust the minute setting in the same way that you adjusted the hour setting. When you have set the time, press the time setting button again as confirmation.

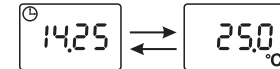3. The display of the device now alternates between the time and the ambient temperature.

Using the device as a room thermometer and clock constitutes a continual measuring process, which significantly reduces the battery service life.

Operating Thermometer

Bear in mind that the thermometer needs to have been in the room in which the measurement is taken for at least 30 minutes before use.

- Some people produce different readings in their left and right ear. In order to record temperature changes, always measure a person's temperature in the same ear.
- The ear thermometer may be used by children only under adult supervision. Measurement is usually possible over the age of 6 months. In infants under 6 months, the ear canal is still very narrow so the temperature of the eardrum often cannot be recorded and the result displayed is often too low.
- The measurement must not be taken in an ear affected by inflammatory diseases (e.g. discharging pus or secretion), after possible ear injuries (e.g. eardrum damage) or in the healing period after operative procedures. In all of these cases, please consult your doctor.
- Use of the thermometer on different persons can be inappropriate in the event of certain acute infectious diseases because of the possible spread of germs despite cleaning and disinfection. If you have any doubts, please consult your doctor.
- This thermometer may only be used without a disposable protective cover.
- If you have been lying on one ear for some time, the temperature is slightly raised. Wait a little while or measure in the other ear.

• As ear wax can affect the measurement, you should clean the ear before measuring if necessary.

Measuring body temperature in the ear

1. Press the Φ/\rightarrow button for 1 second to switch the thermometer on. Following a successful self test, the device emits two short beeps.

2. It is recommended to use your thumb to pull from either side of the Front Cover to remove it easily.

3. Make sure that the sensor tip and also the ear canal are clean. As the ear canal is slightly curved, you have to pull the ear slightly up and backwards before inserting the sensor tip.

4. Insert the sensor tip carefully and press the Φ/\rightarrow button for 1 second.

5. Release the Φ/\rightarrow button. The end of the measuring time is signalled with a short beep and the measured value appears on the display.

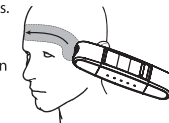
Measuring body temperature on the forehead

1. Please note that the forehead/temple must be free from sweat and cosmetics and that taking medication and skin irritations can distort the result when measuring temperature on the forehead.

2. Press the Φ/\rightarrow button for 1 second to switch the thermometer on. Following a successful self test, the device emits two short beeps.

3. Place the measuring head with the forehead/cover cap fitted on the temples, hold on the Φ/\rightarrow button and move the thermometer smoothly over the forehead to the other temple.

4. Release the Φ/\rightarrow button. The end of the measuring time is signalled with a short beep and the measured value appears on the display.

Battery Installation

When the batteries get weak, the

battery warning symbol appears " H ". It is still possible to measure temperature. When the battery symbol flashes " Lo " appears in the display, the battery must be replaced. If the battery is too flat, then the thermometer will switch off automatically.

NOTE:

- When changing the battery, use battery of the same type, make and capacity.
- Do not use rechargeable battery.
- Use battery free from heavy metals.

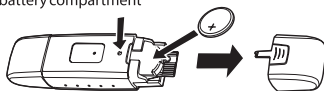
1. Open the battery case by pushing the knob for operating the battery compartment and sliding off the cover on the back of the unit.

2. Remove used battery.

3. Replace with a lithium 3V CR2032 battery in designated area. The larger part of the battery should be facing up. You should still be able to see the "+" sign when battery is installed.

4. Slide the cover back on. Unit is ready for immediate use.

Knob for opening battery compartment

Dispose of used battery in accordance with the applicable legal regulations. Never dispose of battery in the normal household waste.

Trouble-shooting

When a malfunction or incorrect temperature measurement occurs, an error message will appear as described below (comes with backlight and voice pronounces.).

LCD Display:

H Cause: The temperature measured is higher than 43°C (109.4°F)

Lo Cause: The temperature measured is lower than 34°C (93.2°F)

Solution: Operate the thermometer only between the specified temperature ranges. If necessary, clean the sensor tip. In the event of a repeated error message, contact your retailer or Customer Services.

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LCD Display:

 Cause: The operating temperature is not in the range from 16°C~35°C (60.8°F~95°F)

Solution: Operate the thermometer only between the specified temperature ranges

Storage and maintenance

- Clean the probe and tip with an alcohol swab before and after each measurement.
- Use a soft, dry cloth to clean thermometer body. Never use abrasive cleaning agents, thinners or benzene for cleaning. Do not scratch the surface of the probe lens or the display. Do not expose the thermometer to extreme temperatures, humidity, direct sunlight, or shock. Store at room temperature.

Applied Standards

This product conforms to the provisions of the EC directive MDD (93/42/EEC). The following standards apply to design and /or manufacture the products:

- EN12470-5:2003 Clinical thermometers-Part 5: Performance of infrared ear Thermometers (with maximum device)
- ASTM E1965-98 Standard Specification for Infrared Thermometers for Intermittent Determination of Patient Temperature.
- EN980: 2003 Graphical symbols for use in the labeling of medical devices.
- IEC/EN 60601-1 Medical electrical equipment-Part1:General requirements for safety.
- IEC/EN 60601-1-2 Medical electrical equipment-Part2: Collateral standard: Electromagnetic compatibility-Requirements and tests.
- ISO 14971 Medical devices-Application of risk management to medical devices. The Classification according to
- IIEC/EN 60601-1 sub-clauses 5:
 - Internally powered equipment
 - IPX0
 - Equipment not suitable for use in the presence of a flammable anaesthetic mixture with air or with oxygen or nitrous oxide
 - Continuous operation

Technical Specifications

- Range of Measurement : 34°C~43°C (93.2°F~109.4°F)

- Measuring accuracy:

Ear:
±0.2°C(0.4°F): 35.5°C~42°C (95.9°F~107.6°F)
±0.3°C(0.5°F): outside the range
Forehead:
±0.3°C(0.5°F): 34°C~40°C (93.2°F~104°F)

- Operating environment:
16°C~35°C (60.8°F~95°F) with relative humidity up to 85% (non condensing)
- Storage/Transportation environment:
-25°C~-+55°C (-13°F~-+131°F) with relative humidity up to 85% (non condensing)
- Display resolution: 0.1°C or °F
- Power supply: A lithium 3V CR2032 battery
- Weight: approx. 45 g including battery
- Dimensions: approx. 109 mmX30 mmX22mm

Repair and maintenance

Any extensive maintenance or repair that requires access to the inner parts of the product must be performed by an expert service centre.

Environmental protection:

- Packaging materials and obsolete appliances should be recycled.
- The transport box may be disposed of as sorted waste.
- Polyethylene bags shall be handed over for recycling.

Appliance recycling at the end of its service life:

A symbol on the product or its packaging indicates that this product should not go into household waste. It must be taken to the collection point of an electric and equipment recycling facility. By making sure this product is disposed of properly, you will help prevent the negative effects on the environment and human health that would otherwise result from inappropriate disposal of this product. You can learn more about recycling this product from your local authorities, a household waste disposal service or in the shop where you bought this product.

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