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Mi, **Momert Zrt.**

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Operating manual MORE FUNCTIONS DIGITAL SCALE MODEL 5863

Thank you for choosing more functions digital scale. These scales have been designed to accurately indicate your body weight, body fat, body water, muscle, bone weight, BMI (Body Mass Index) and daily calorie need. These scales should give many years of service under normal use.

Please study the entire Operating Manual carefully before you start using the product. Keep the manual in a safe place for future reference. Make sure other people using the product are familiar with these instructions.

SPECIFICATION:

- Precision measuring system
- Automatic zero resetting
- Automatic switch-off
- Low power indicator
- Overload indicator

Technical parameters	
Voltage	3 pcs 1,5 V AAA batteries (included)
Weighing range	10-180 kg
Weighing increments	100 g
Weighing units	kg/lb/st
Display size	LCD diameter= 10 cm/ 3,9 inch
Body fat division	0,1%
Water division	0,1%
Age range for body fat	10-100 years
Height range	100-225 cm
Height unit	cm/ft/in
Memory	10 person

RECOMMENDED WORKING ENVIRONMENT:

- Temperature: 0 – 40 Celsius / 32- 104 Fahrenheit
- Relative humidity: < 85%

IMPORTANT SAFETY INSTRUCTIONS:

- Make sure you follow the marked polarity of battery when inserting and replacing it.
- If the appliance will be out of use for a longer time, remove the battery.
- Do not leave the unit unattended when in use.
- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
- Use the product only for the intended purpose specified in the Operating Manual.
- Use the unit only on a stable horizontal surface.
- Before cleaning, assembly, or disassembly of the cover, switch the unit off and remove the battery!
- Do not use abrasive or chemically aggressive substances for cleaning the unit. Clean after use with a lightly damp cloth, do not use solvents or immerse the product in water.
- Do not use the unit if it does not operate properly, or if it has been dropped, damaged or exposed to any liquids. Take the unit to an authorized service center for testing and repair.
- The unit is suitable for domestic use only. It is not designed for commercial operation.
- Keep the scales in a dry place to protect the electronic components.
- Keep the unit away from heat sources, such as radiators, ovens etc. Do not expose the unit to direct sunlight or moisture.
- Do not touch the unit with wet or damp hands.
- The surface of the scales will be slippery if wet.
- Weigh yourself without clothing and footwear, before meals and always at the same time of day. For better accuracy, it is not advised to take the measurement within 2 hours of waking up.
- The product is a fragile measuring device. Prevent it from crashing or falling.
- Avoid overloading the scales, otherwise, permanent damage may occur.
- Do not repair the appliance yourself. Contact an authorized service center. Do not disassemble the scales other than to replace the battery.

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
- A user who has any metal material inserted into his/her body may get differences in measurement of body fat reading.
- Pacemaker users are advised against using this device.

FAILURE TO FOLLOW THE MANUFACTURER'S INSTRUCTIONS MAY LEAD TO REFUSAL OF WARRANTY REPAIR!

INSTRUCTION MANUAL:

1. The scale includes replaceable lithium battery. Please remove the protective film of the battery before use.
2. Place the scale on a hard, flat surface. (Do not use it on carpet or soft surface)
3. Tap on the top platform of the scale to power on. Wait a few seconds until the following will be displayed: 0.0
4. Select the weighing mode using the switch (kg/st/lb) located on the back of the scales.
5. Step onto the scale carefully. Your weight should be equally distributed on the scale. Stand still, while the scale measures your weight. The result, your accurate body weight will be displayed in the following way flashing twice on the display.
6. When you step off the scales it will automatically switch off after 8 seconds.
7. If no weighing is carried out while 0.0 is displayed the scales will automatically switch off after 8 seconds.

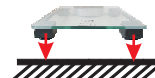BODY FAT MEASUREMENT:

In order to ensure accurate results users should enter the required personal parameters: age, sex and height.

These scales have the facility to store the personal parameters of 10 people.

The technology used in these scales, Bio-electrical Impedance Analysis (BIA) is more technologically advanced than the older Body Mass Index (BMI) method where fat % is calculated using only height and weight. BIA technology uses 5 parameters to get a more accurate resolution for the analysis

(weight and biological resistance which are automatically measured by the scales together with the three personal parameters entered manually).

MODE SELECTION AND PARAMETER SETTING:

- Press and hold **SET** to switch on the scale.
- Press **▲** or **▼** to switch between body parameter 0-9. Press **SET** to confirm.
- Press **▲** or **▼** to select male or female. Press to confirm.
- Press **▲** or **▼** to indicate age. Press **SET** to confirm.
- Press **▲** or **▼** to indicate height. Press **SET** to confirm.
- Select the weighing mode using the switch (cm/ft-inch) located on the back of the scales. Press **USER** to confirm.

The values are flashed and recorded. Stand on the scales after 0.0 appears. First the weight is recorded, then SSS is displayed, indicating that the scales are analysing, please do not alight. Once finished, the weight, body fat , body water, muscle, bone, BMI and KCAL will be displayed alternately for 2 times before the scales automatically turn off. At the next measuring, if you turn on the scale with SET button, then choose your user ID and push USER button, then the scale shows the previously measured datas. If you only want to measure your weight, then tap on the scale. Compare the body fat reading against the (Body fat and Body wter ratio Chart) to confirm the body status.

BODY FAT RATIO CHART %				
FEMALE				
AGE	Low	Healthy	High	Obese
< 39	< 21	21-33	33-39	>39
40-59	< 23	23-35	35-40	>40
>59	< 24	24-36	36-42	>42
MALE				
AGE	Low	Healthy	High	Obese
< 39	< 8	8-19	19-25	>25
40-59	< 11	11-22	22-28	>28
>59	< 13	13-25	25-30	>30

WATER & MUSCLE RATIO CHART %								
AGE	FEMALE				MALE			
	Low	Healthy	High	Muscle	Low	Healthy	High	Muscle
10-15	<57	57-67	>67	>34	<58	58-72	>72	>40
16-30	<47	47-57	>57		<53	53-67	>67	
31-60	<42	42-52	>52		<47	47-61	>62	
61-80	<37	37-47	>47		<42	42-56	>56	

BMI (BODY MASS INDEX)								
FEMALE				MALE				
Low	Healthy	High	Obese	Low	Healthy	High	Obese	
<18,5	18,6-24	24,1-28,9	>29	<18,5	18,6-25	25,1-29,9	>30	

BONE WEIGHT CHART (Kg)						
	FEMALE			MALE		
	Weight	< 45 kg	45~60 kg	> 60 kg	< 60 kg	60~75 kg
Bone weight	~ 2,6 kg	~ 3,2 kg	~ 8 kg	~ 3,4 kg	~ 3,8 kg	~ 8 kg

KCAL						
AGE	FEMALE			MALE		
	Weight (Kg)	KCAL	Weight (Kg)	KCAL		
15-17	50	2116	60	2525		
18-29	55	2136	65	2502		
30-49	60	2094	70	2489		
50-69	60	1948	70	2278		
70-80	50	1723	65	1956		

ERROR WARNING:

In case of error (e.g.: overloading), the following is displayed:

**CHANGE BATTERY:**

The battery needs to be replaced, when the display screen shows:

Please open the battery holder at the backside of the scale and carefully take out the battery by pressing the small silver colour holding plate to get the battery out of the battery holder.

Replace the low battery to a new battery. The low battery should be thrown to the appropriate place.

BATTERIES ARE NOT INCLUDED IN WARRANTY!**CLEANING AND MAINTENANCE:**

The surface of the unit may only be cleaned using a wet cloth. Do not use any detergents or hard objects as they may damage the unit's surface!

Never clean the unit under running water, do not rinse it, and do not submerge it in water!

REPAIR AND MAINTENANCE:

Any extensive maintenance or repair that requires access to the inner parts of the product must be performed by an expert service centre.

ENVIRONMENTAL PROTECTION:

- Packaging materials and obsolete appliances should be recycled.
- The transport box may be disposed of as sorted waste.
- Polyethylene bags shall be handed over for recycling.

Appliance recycling at the end of its service life:

A symbol on the product or its packaging indicates that this product should not go into household waste. It must be taken to the collection point of an electric and equipment recycling facility. By making sure this product is disposed of properly, you will help prevent the negative effects on the environment and human health that would otherwise result from inappropriate disposal of this product. You can learn more about recycling this product from your local authorities, a household waste disposal service or in the shop where you bought this product.

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